



## Install the SmartMan Program

This program will run on any Windows or Vista PC. It will not run on under the Macintosh O/S.

The file is an executable self-extracting file that will:

- 1) Create the v4ems subdirectory in the Program Files area on your C: Drive (c:\v4EMS).
- 2) Install the program into this directory.
- 3) Create the data directory (c:\v4EMS\_DATA). Every login that is created will place a record in this directory.
- 4) Put an icon on your desk top.

To install the latest Version you can click on this [UPDATE LINK](#) and install from there.

As an alternative, install from the CD and click on Upgrade Software from the Main Menu.

You must be connected to the Internet to use either of these options.

You can also use this link to [Download latest version of SmartMan](#). Here is the direct link address: ([http://www.v4ems.com/Mods/V4EMS\\_Update.html](http://www.v4ems.com/Mods/V4EMS_Update.html)).

If you are upgrading and get an error when trying to install, go to the c:\program file\v4EMS subdirectory and delete it. Try to install the program again.

## Starting the Program



Once you have installed the program click on one of the SmartMan icons (on the desktop or in the start menu) to start the program.

1. **Adapt the Screen or Not:** The first screen gives you the option to adapt the size of the screen or to continue at the current resolution. If “Adapt” is selected the screen will return to the original resolution when the program is closed but you must make sure that you exit the program from the ‘Quit’ button at the top of the screen then click ‘Yes’ button in the ‘Are you sure you want to quit’ box.
2. **Login:** You must create a login name. This can be a short name or a nick name. You will use this name each time you go into the program.
3. **Password:** Because every person’s scores are private, when a person registers he or she must give their login name a password. This must be at least six letters or numbers long and cannot have any spaces. When you have typed your password, click continue. You will be asked to type your password a second time. NOTE: both the ‘Log In’ and ‘Password’ are case sensitive.
4. **Connection Verification:** SmartMan will now check the connection to the manikin.
  - a. **Good Connection** - If the connection via all cables is okay, you will see a popup window. If this is the first time you logged in to SmartMan, it will explain that you MUST register. Click “Continue” to be taken to the main menu.

- b. **Bad Connection** – If the connection cannot be established, a popup window will explain that an error has occurred. The connection status will indicate “-1”. Check that all cables are correctly plugged in. Then click the “Restart” Button. This will take you to the adapt screen page and you will proceed through the beginning steps again.
1. **Main Menu:** There are three menus on the Main Screen. All buttons must be white in order to be active. In the example below, the “register” Button is not available because the person has already registered. All of the other options are available.

### SKILLS MENU


This is where you select which skills to practice. When you click on one of the selections a drop down menu will give you further options.

### ACTIONS MENU

This is for registering, changing the person who is logged in, looking at previous results, connecting to the SmartMan Home Page, Checking for Updates and Calibrating SmartMan

### CONTENT MENU

This explains to a new user how to perform some of the basic skills

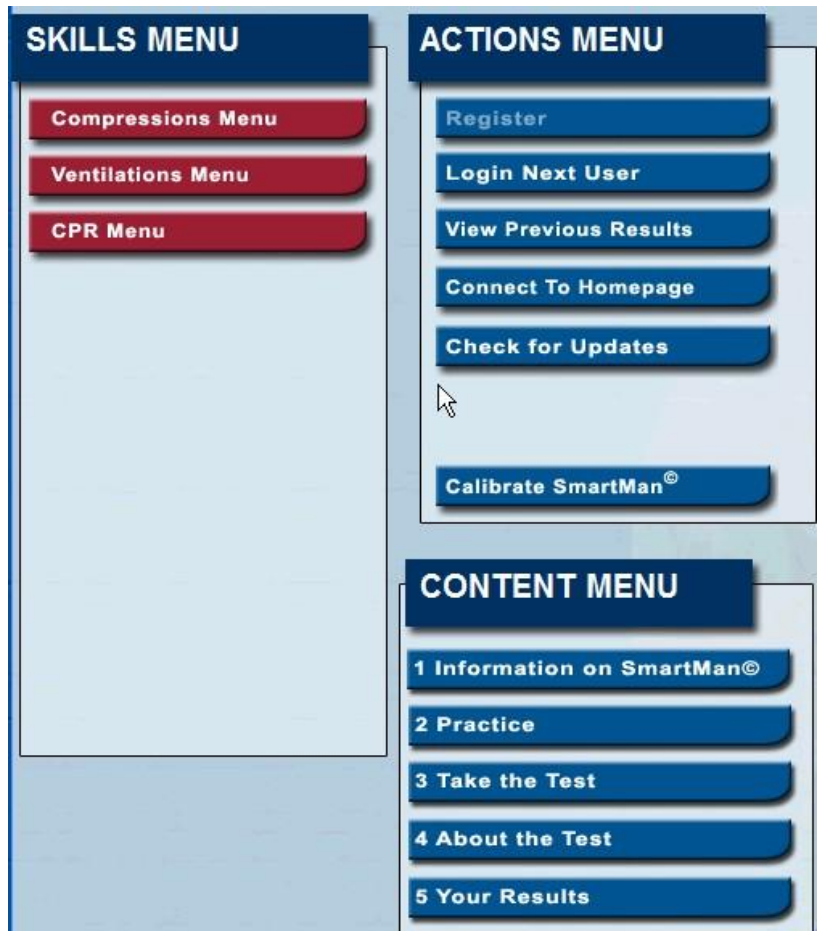

2. **First Time Log In:** If this is the first time you logged in, you must **register**. Click on the “Register” Button at the top of the Actions Menu. For Corporate set up, Corporate registration and Individual registration instructions please proceed to the ‘Registering’ section of the instruction pages.